





Scope-Shield and Scope-Cone

The Nuview Scope-Shield and Scope-Cone have been developed in conjunction with leading dental clinicians to offer shielding against spray and droplets during routine microscopically assisted dental procedures.

Scope-Shield

The Nuview Scope-Shield is a unique two part system which ensures the closure of gaps at the binocular/shield interface. The two-part system with the rotatable baffle caters for the majority of inter-pupillary distances. The shield ocular diameter is designed for microscopes with a minor ocular diameter of 43mm. If an assistant scope is used a further Scope-Shield can also be utilised.

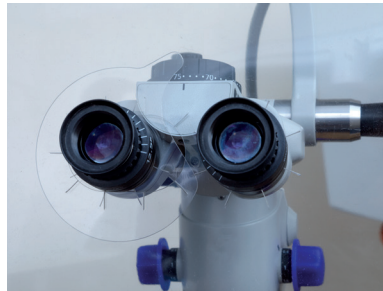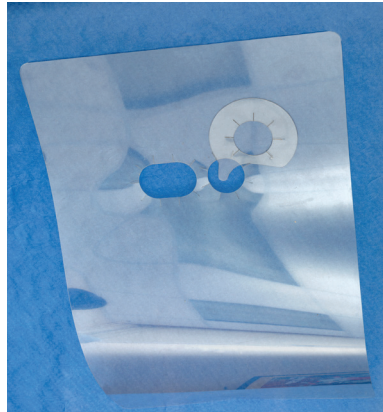


Scope-Cone

The Nuview Scope-Cone shields the microscope head from spray generated during dental procedures. The conical nature of the shield ensures access to the microscope handles during use. A lens cap system is used to secure the Cone to the microscope.

Both the Scope-Shield and Scope-Cone are designed to be reusable but it is recommended that they are replaced on a daily basis. Both products are provided in packs of 5. Cleaning and disinfecting should be undertaken in compliance with the practice/institution approved infection control procedures.

The Nuview Scope-Shield and Scope-Cone are intended to Guard the microscope and user against droplets and spray - they are not intended to be used as a substitute for personal protective equipment.

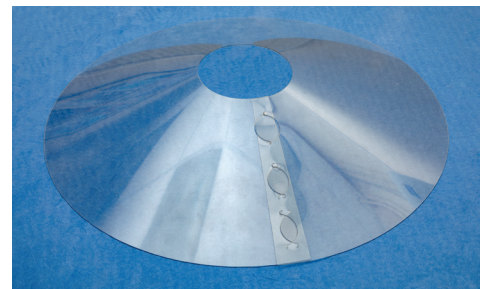

Order details

Scope-Shield & Baffle (Box 5)	PPS-1001
Scope-Cone Lens Cover (Box 5)	PPS-2002